

Precise dissection in thin-walled anatomy

Benefit from an all-in-one solution uniting cutting-edge technology with new flexibility.

Inject, cut, coagulate – accelerate your workflow by instant application without changing instrument.

TAILOR YOUR OWN TECHNIQUE.





Elevating 3rd space endoscopy

Fine electrode for precise resection

Coaxial nozzle for high-pressure injection

Flexible probe for smooth passage even in retroflexion

Small diameter for adequate suction through the working channel

Locking mechanism for constant electrode length



BENEFIT FROM AN ALL-IN-ONE SOLUTION WITH HYBRIDKNIFE® FLEX AND THE GI WORKSTATION

- ☑ Needle-free injection through mucosa
- ☑ Repeated Micro-cushioning
- ☑ VIO® modes (e.g., preciseSECT®)
- ☑ Minimal thermal spread
- ☑ Effective coagulation and hemostasis

Inject needle-free


Expand the submucosal space as needed and without changing your instrument using the high-pressure, needle-free, tissue-selective hydrodissection.

Cut precisely

The thin and short electrode enables high electrosurgical precision with minimal thermal damage.

Coagulate as needed

Apply various coagulation techniques by adjusting tissue contact and using VIO® electrosurgical modes

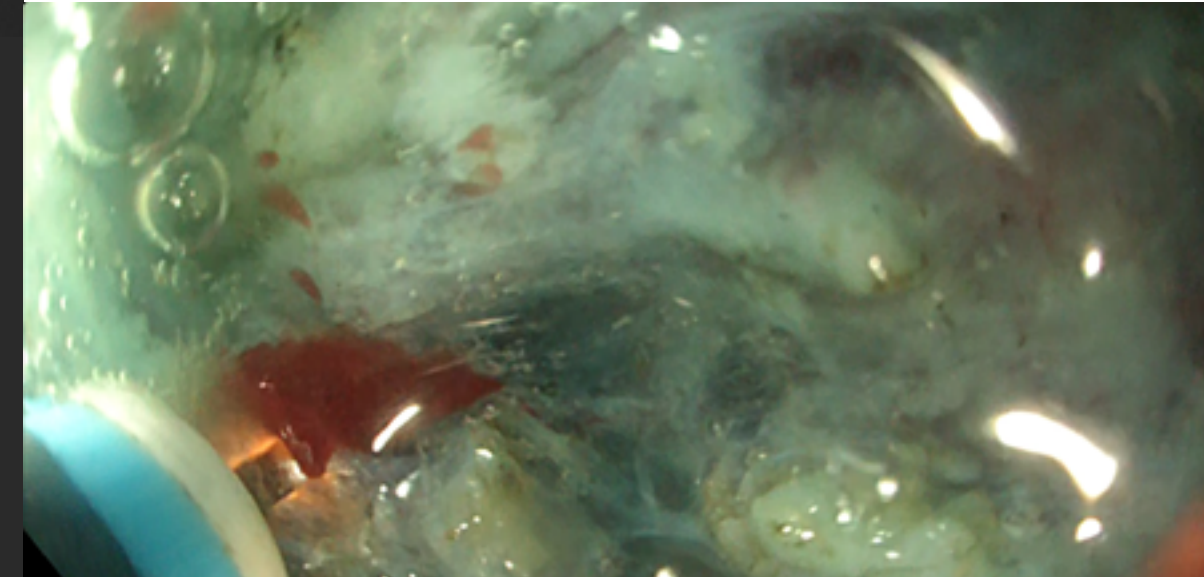

Tubing length	2300 mm
Outer diameter (distal end)	2.6 mm
Outer diameter (tubing)	2.1 mm
Electrode length	1.5 mm / 2.0 mm
Type	I
Retro-flexibility	More than 210°
I-Type 1.5 mm	No. 20150-112
I-Type 2.0 mm	No. 20150-113

Marketing claim list 0188716, V005